

Flag-Tag Mouse Monoclonal Antibody

Catalog Number: AT013

Product Package

Cat. No.	Unit Size	Concentration	Storage
AT013	100 μ L	1 mg/mL	Short term: 4°C Long term: -20°C Avoid freeze/thaw cycle

Introduction

Flag-Tag antibody is part of the Tag series of antibodies, the best quality in the research.

Storage Condition and Buffer: 1 mg/ml in PBS, pH 7.4, containing 0.02% sodium azide and 50% glycerol. Store at +4°C short term. Store at -20°C long term. Avoid freeze / thaw cycle.

Application: WB ICC/IF IP.


Reactivity: All

Dilution: WB 1:3000~1:10000 ICC/IF 1:100~1:200 IP 1:100.

Purification: Affinity-chromatography.

Specificity: Flag-Tag Mouse Monoclonal antibody detects endogenous levels of C-terminal, internal, and N-terminal Flag-tagged proteins.

Immunogen: A synthesized peptide.

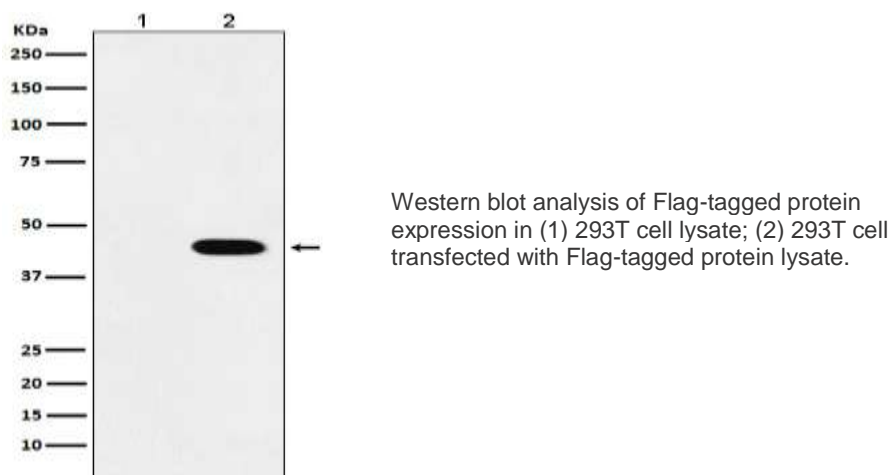

For Research Use Only.
www.abpbio.com